


Dilatometer sample holder for powders

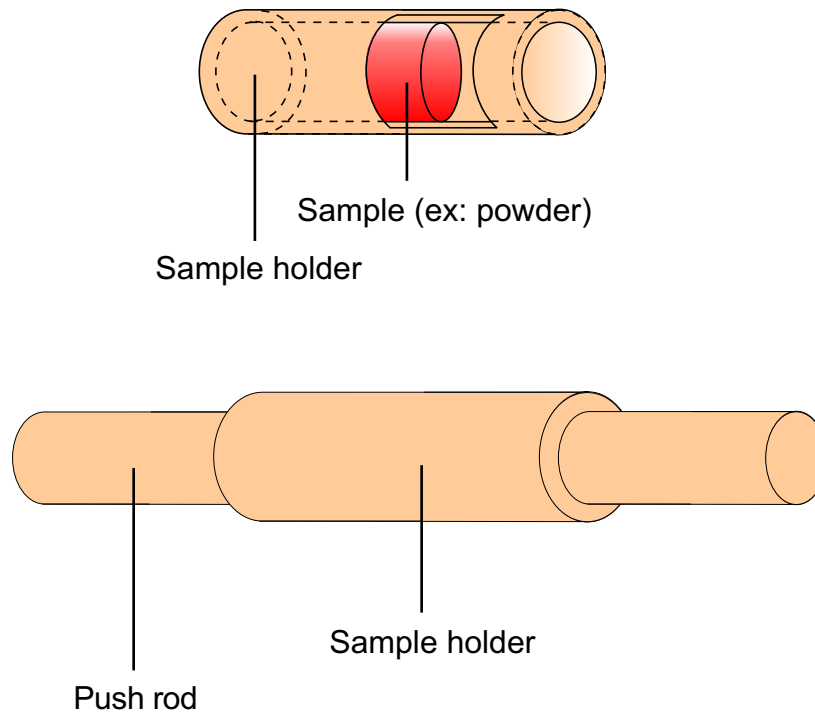

In response to customer requests to test powder samples in the dilatometer, Linseis has developed a new powder sample holder. Using the new powder sample holder it is now possible to test powders and examine binder burnout, CTE, and sintering using a standard sample holder along with the optional powder sampler holder.

The powder sample holders are available made from Al_2O_3 , quartz, and graphite.

There are two versions available, one with fixed end and movable push rod. The other version has two movable push rods inserted in the holding sleeve.

The powder container can be used in horizontal or vertical dilatometers.